

Implantech[®]
Superior Patient Aesthetics

NEWS FOR YOUR PRACTICE

AUGUST 2022
800.733.0833 | implantech.com

SPOTLIGHT ON Vertical Mandibular Angle

For lowering the jaw angle

PICTURES SPEAK VOLUMES

Check out some before-after photos

TRAINING WORKSHOPS

Sign up for updated class schedules

FANCY MEETING YOU

Meetings where we hope to see you!



WELCOME TO Fabulous LAS VEGAS NEVADA

MASTERS EDUCATIONAL Series COURSES

Facial & Body Contouring Implants

Saturday & Sunday
NOVEMBER 12-13, 2022 **Implantech**

Advanced Techniques Using Facial & Body Contouring Implants

When's the last time a weekend added a new dimension to your practice?

Give your patients a great set of options by attending our two-day **Masters Educational Series** workshop in Las Vegas.


The Vertical Mandibular Angle

Literally, a jaw-dropping look!

Traditional mandibular implants were designed to sit on the patient's jaw angle to increase its lateral projection. This can be aesthetically effective, unless the jaw angle is within 2 cm of the earlobe.

An implant on such a high jaw angle can make the patient's face look too full or chunky, instead of creating a stronger, more defined jawline.

The Solution – Implantech's **Vertical Mandibular Angle (VMA)**.

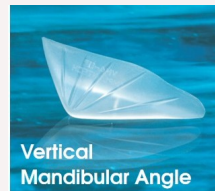

This innovative silicone implant is designed to lengthen the ramus (i.e., increase the vertical distance of the jaw angle from the earlobe).

It then creates a distinct, more angular jaw shape below the patient's natural jawline.



The Placement –

The VMA is placed in a subperiosteal pocket that should follow the outline of the mandibular ramus (vertically), the lateral

ramus posteriorly, and the lower border inferiorly.

Screw fixation is particularly important where a portion of the implant is not on the bone, as the VMA must resist forces of displacement caused by masseter muscle

With clear, thorough instruction, our expert faculty imparts the latest techniques using facial and body-contouring implants made of supple, solid silicone.

Your Esteemed Faculty



William J. Binder
MD, FACS



Barry L. Eppley
MD, DMD



Stanley A. Okoro
MD, FACS




Rady Rahban
MD, FACS



Robert J. Troell
MD, FACS


You'll learn how to achieve aesthetic and long-lasting patient outcomes from both classroom and cadaver lab sessions.

Your topics covered each day include:



FACIAL
DAY 1 **HANDS-ON**
CADAVER LAB

Implants covered
Chin • Jaw Angle • Malar • Submalar
Paranasal • Nasal • Temporal Shell

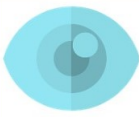

BODY
DAY 2 **OBSERVE**
CADAVER LAB

Implants covered
Gluteal • Pectoral • Calf

Also discussed on both days:

- Patient-specific implants
- Challenging cases
- Managing complications



Upon course conclusion, you'll receive a Certificate of Completion, plus a "Course Certified" icon for prominent display on your website.

contraction.

The Result (for him and her!) –



Contact Implantech for more information about the Vertical Mandibular Angle implant or any of our superior products.



Did you
KNOW?

You can place your orders online!
Just register on our website (USA only)

Questions
about
our products or
their uses?

800.733.0833

Implantech.com

Attend one or both days, but make sure you take advantage of this opportunity!

RESOURCES



USEFUL INFORMATION FOR
YOU AND YOUR PRACTICE

IMPLANTECH | implantech.com



IMP-642/0822